



Chromai[®]

— 科诺美 —



Leaps 2D

Two-Dimensional
Ultra-High Performance
Liquid Chromatography System
UHPLC



Innovation for Better Life

Leaps 2D

Two-Dimensional Ultra-High
Performance Liquid
Chromatography System

Stay one step ahead in
complex analysis


Flexible
Configuration

Automated
Operation

High Column
Capacity



Precisely meets customers' differentiated
analytical needs

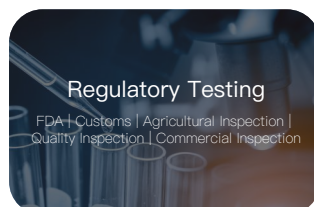

A pathway to unlock complex application scenarios

Simplify pretreatment and multiply analytical efficiency


- Advanced UHPLC analytical technology significantly boosts analytical efficiency
- Online cleanup and enrichment simplify pretreatment and greatly improve detection sensitivity
- Meets the needs of analyzing complex matrices such as food matrices, natural products, and therapeutic drug monitoring samples

Multi-column switching for ideal measurement of complex samples

- With a multi-column switching system, low-concentration samples can be measured and multiple sample types can be analyzed in parallel
- Supports automated transfer between first-dimension and second-dimension columns through automatic switching
- Fully automated sample loading greatly improves usability and eliminates complicated manual intervention
- Analytical results can be rapidly obtained via heart-cutting, multiple heart-cutting, and high-resolution sampling
- Easily enables automatic selection between 1D and 2D modes, choosing the appropriate separation performance based on sample complexity

Flexible options with perfect upward and downward compatibility

- Interference-free method transfer, seamless upgrades, and worry-free operation
- Excellent module compatibility allows easy upgrades and protects your investment
- Parallel dual columns and series SPE greatly improve efficiency
- Flexible and powerful dual ternary pumps support parallel and independent operation

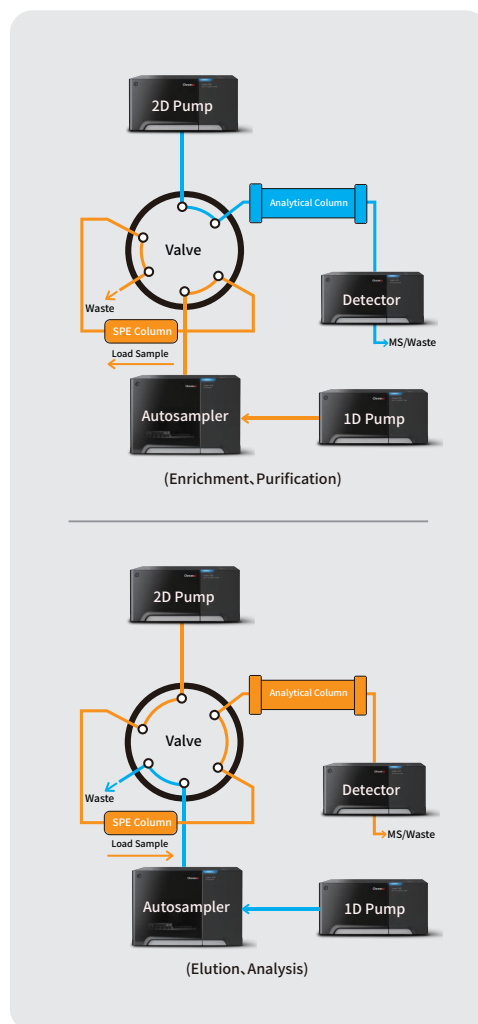


Superior Performance, Reliable Results

● Diode Array Detector (DAD)



- Pixels: 1024 pixels
- Wavelength Range: 190–800 nm
- Slit Width: 1.2 nm / 6 nm
- Monitoring Channels: 12 channels simultaneously
- Additional Functions: 3D spectra, peak purity
- Flow Cell: Max path length 60 mm

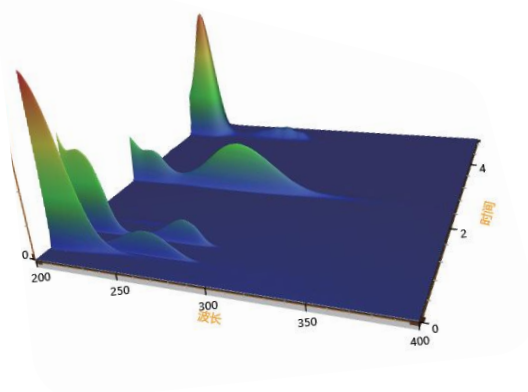


● UV-Vis Detector



- Optical Design: Automatic switching
- Wavelength Range: 190–900 nm
- Acquisition Rate: Up to 100 Hz
- Linear Range: 2.5 AU
- Standard Features: Wavelength time program, switching, stop-flow scanning
- Maintenance: Easy lamp replacement

● Fluorescence Detector

- Sensitivity: New optical design
- Stray Light Suppression: Excellent S/N
- Wavelength Range: 200–750 nm / 900 nm
- Sampling Rate: Up to 100 Hz

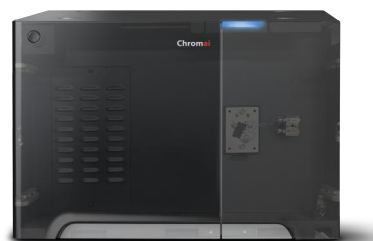

EyouLab Intelligent Data System

Enhances laboratory efficiency with greater accuracy and speed
Effectively reduces operational costs and risks



EyouLab CDS

EyouLab fully complies with FDA 21 CFR Part 11, GMP, GLP, and other regulatory requirements. It offers multiple deployment options, including standalone, networked, and medical versions, with multi-lingual support covering everything from sample preparation to precision medicine. Built on dozens of proprietary technologies, it combines a “smart guided interface + cloud-based data management” to simplify complex experimental operations, ensure data integrity and accuracy, and significantly boost lab capabilities and efficiency—accelerating the transition to intelligent laboratories.




Control Panel

Multi-channel Data Acquisition

Automated Sequences

Instrument Monitoring

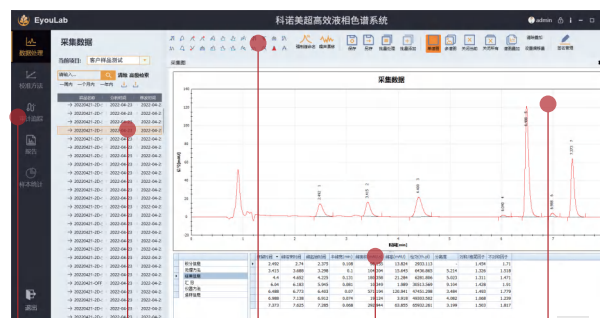

Data File List

Complete Attribute Parameters

Powerful Data Processing

Multi-mode Data Display

Standalone Operation

Remote Access

Archiving Management

Audit Trail

Risk Management

Unified Data Management and Retrieval

Full Regulatory Compliance

Ensures Physical Security of Electronic Data

Easy Remote Access Control

Open-access Network Analysis System

Simple Installation, Startup, and Maintenance

Each acquisition terminal can control up to 4 systems



Leaps

UHPLC System



Leaps Pro

High Throughput UHPLC System



Frontier

UHPLC System

Analytical



Leaps 2D

2D UHPLC System

2D System



Leaps Pro Bio

Bio-Inert LC System

Bio-Inert



Leaps Prep

Preparative LC System

Preparative

Website: www.chromai.com

Chromai · Beijing

Address: 5th Floor, Building 10, Liandong U Valley · Lanbei
Technology Park, No. 3 Xixing Road, Shunyi District, Beijing, China

Tel: (86-10)80429960

Chromai · Suzhou

Address: 5th Floor, Building 5, No. 1268 Jiaotong Road, Wujiang
District, Suzhou, Jiangsu, China

Tel: (86-512)63226066

Email: contact@chromai.com



WeChat
Official Account



LinkedIn
Official Account